

Eleganza 1

официальный дистрибьютор LINET
в Российской Федерации
ООО "МЕДЭК Инструмент Групп"
119602, Москва, ул.Коштоянца 20/4
тел.: +7(495)755-55-31, 755-55-32
факс +7(495)755-55-30
mail@linet-ru.com, medec@medec.ru
WWW.LINET-RU.COM



excellent solution for a good price

Eleganza 1





- Ergoframe® against pressure ulcers
- Easy mobilisation
- Available with a number of controls



Patient

- Easy clean concept
- Safety features such as the Extender or Protector
- Easy bed manoeuvrability



Staff

- Low TCO – total cost of ownership
- Universal bed for a low price
- Overload protection



Management

Improved patient care

COMFORT FOR EVERYONE

The bed accommodates hospitalised overweight or taller patients.

- Robust construction.
- Safe working load 180–250 kg.
- Mattress platform width 90 cm.
- Possible bed extension by 15 cm.



CONTROLS WITHIN REACH

Several clear and easy to use controls are immediately at the patient's disposal.

- 1 Handset**
to control the bed positions can also come with a back lit keyboard.
- 2 Satellite control panel**
comes on a flexible arm and is easy to access, the buttons are embossed for vision impaired patients.
- 3 Controls in the siderails**
are easy to use and always to hand.





ERGOFRAME® +16 cm

When repositioning the backrest and thigh rest, the Ergoframe® function minimises the pressure and shear forces in the patient's pelvic region. This results in a significant reduction in the development of pressure ulcers.



EASIER MOBILISATION

The mattress platform is equipped with an integrated holder for the MobiLift® which supports independent patient mobilisation.

- Respects the different height of patients and can be set at two height levels.
- Always set at the optimum supporting exit point.



Safety first

SAFE SIDERAILS

The bed can be equipped with plastic split siderails or single $\frac{3}{4}$ siderails. In each case they comply with valid standards and all gaps on them are designed with regard to safety.

- The safety of the collapsible plastic siderails is guaranteed by a system of two independent unlocking movements.
- The single siderail controls are located ergonomically for the first time at the top of the siderail.



THE UNIQUE PROTECTOR

The use of the Protector enables protection along the entire length of the mattress platform. The Protector is located at the foot end of the bed and is linked to the single collapsible siderail both in terms of dimensions and design.





BETTER PROTECTION

In some cases it is appropriate to provide the patient with a higher level of protection. An easy-to-assemble extension piece makes the siderail 10 cm higher. The Extender is especially recommended when using an active mattress – up to 26.5 cm height.


GO BUTTON

The green GO button protects the patient against an injury caused by any unintentional movement of the bed. Use the GO button to activate the bed control units for a period of 3 minutes from when last pushed.



For easy hygiene



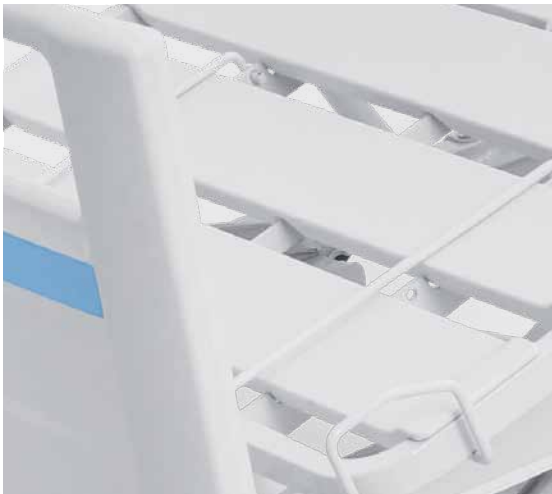

PLASTIC MATTRESS PLATFORM

The mattress platform consists of eight separate removable plastic parts and is produced to make bed hygiene very easy and fast.

- Each part weighs only 700 g, which facilitates hygiene.
- Embossed patient pictogram to identify the part's correct position.



JOINT-FREE TECHNOLOGY

The latest technology is used during the manufacture of individual parts to minimise the space where unwanted infections may hide. In terms of easy cleaning, the bed is very safe and user friendly.



SIMPLE CONSTRUCTION

From the start to achieve easy cleaning, the Eleganza 1 bed was designed to maximise the simplification of the structure to include individual parts and surfaces.

Staff friendly





5th CASTOR

The activated 5th castor makes transfer through long corridors easier. Creating a turning point in the middle of the bed, it allows for tight turns, also providing significant reduction in the force needed to turn the bed.



SUPERVISOR CONTROL PANEL

The panel has a coiled cable providing flexibility of movement and operation of the bed from an optimum position.



CENTRAL CONTROL PANEL

The panel for the control and deactivation of the functions is located under the foot end and is integrated into the bed frame.




FOOT CONTROLS

The controls allow the practical and easy adjustment of the height and examination position. They are equipped with a mechanical protector against accidental activation.

Eleganza 1 Extra SAFE

CARE AND ABSOLUTE SAFETY

The Eleganza 1 Extra SAFE offers a high quality standard of hospital bed. Thanks to the modifications it is suitable for specific departments such as psychiatric, prison hospital, mental health etc.

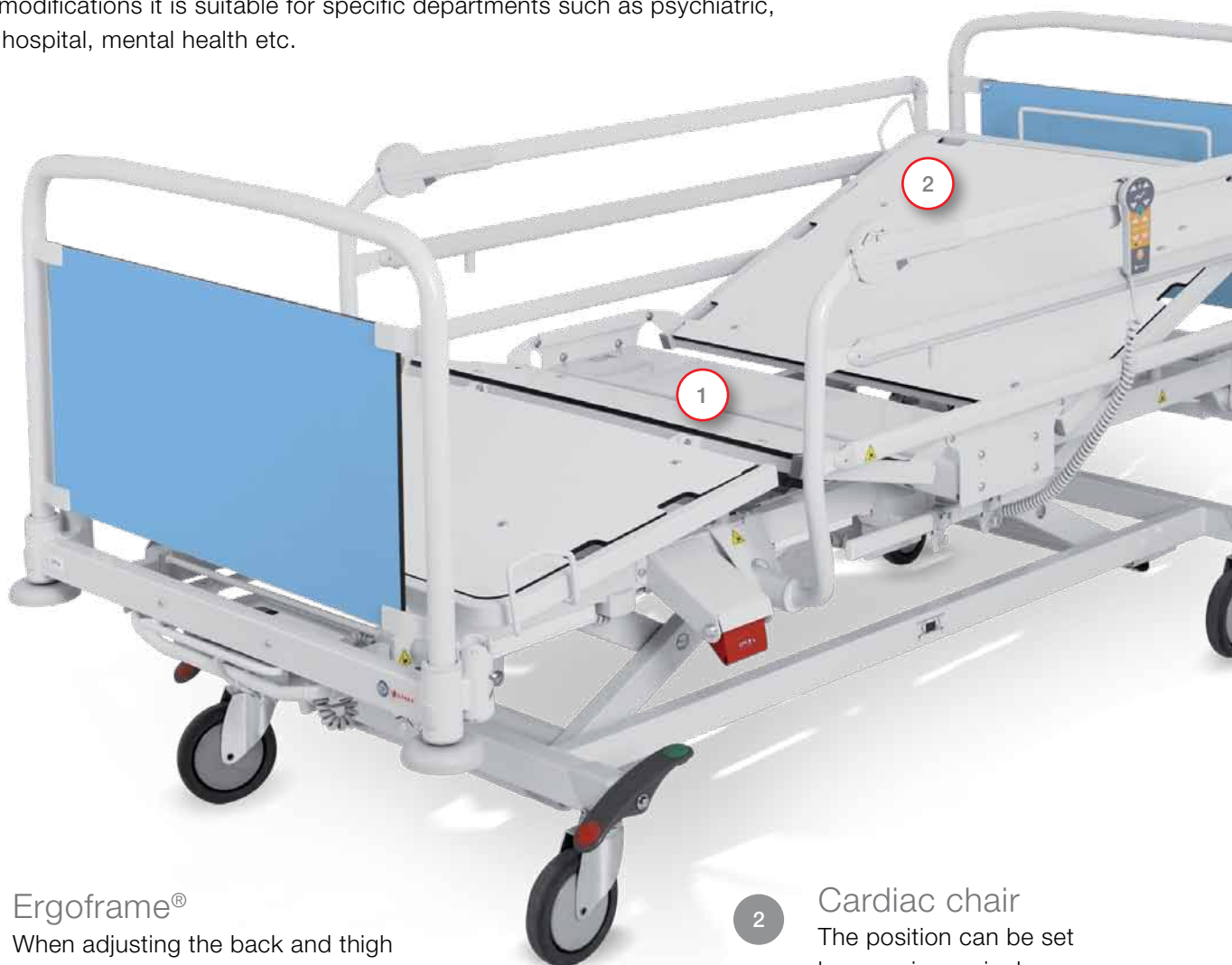

1 Ergoframe®
When adjusting the back and thigh rest the Ergoframe® system extends the mattress platform by 16 cm.

2 Cardiac chair
The position can be set by pressing a single button.

Basic equipment of Eleganza 1 Extra SAFE *

Fixed mattress platform covers	Plug and Play connector of the patient's handset
Fixed tube ends	Back-up battery
Fixed end decors	Supervisory panel with cardiac chair functions, TR and ATR by pressing a single button
Built-in supervisory panel in the line shelf	Protection of the supply cable in a locking box

* Unlike the standard Eleganza 1 version



FIXED PARTS

The entire bed and its panels are firmly installed and the mattress platform parts are also fixed. All smaller parts of the bed, such as the hooks on the accessory bar, are secured.



CABLE PROTECTION

The power supply cable is in a box with a lock. The cables are fed in a way so they do not endanger the patient.



SEGUFIX KEY

The Segufix is compatible, and the operator does not need any additional keys.



FIXED VERSION

A bed in a version without castors can be fixed to the floor.

Eleganza 1 Nursing



In its nursing-care version the Eleganza 1 increases the quality of long-term care. The function and properties fully comply with the needs of nursing homes and long-term and after-care facilities.

**Excellent
DESIGN**

WOODEN DECORS

The colour scheme of the bed is highly variable and fully conforms to the environment of the specific facility. By choosing his own design and version the customer will acquire a functional and totally individual model.



R 5184 Royal Maple



R 5320 White Beech



R 5313 Planked Beech



R 4601 Golden red Alder



R 5360 Precious Cherry



R 4801 Standard Walnut



R 5681 Havanna Cherry



Eleganza 1 with full-length siderails.

Eleganza 1 LE

COST-EFFECTIVE VERSION

The Eleganza 1 LE bed construction is the same as the classic Eleganza 1 model: lift, maximum load, undercarriage including all versions of castors and brakes, mattress platform and ends. It also comes with either single collapsible siderails or split plastic siderails.

The bed's control unit is simpler.





The simplified version of the Eleganza 1 LE bed offers high functionality at a more affordable price.

- Universal control box 100–240 V.
- Bed can be equipped with a Plug and Play handset.
- 1 lever of mechanical CPR (possible upgrade for standard 2 levers).



BED CONTROL

The Eleganza 1 LE is equipped with a simpler control unit, so the bed can be controlled by a single handset. A simple handset controls the electrical positioning and height setting.

- With mechanical lock.
- 8 or 10 buttons (two buttons for tilting into the TR/ATR position).



Eleganza 1 LE nursing version.



A more affordable version of the non-tilting electric bed.

Practical details

1 Removable bed ends can be fixed to the mattress platform frame with simple eccentric locks.

2 Horizontal corner bumpers protect the bed.

3 Mattress platform is easily removable and cleanable and is compatible with 86 or 90 cm mattresses.

4 Linen shelf may also be used to store the supervisor panel.

5 Mechanical CPR – additional option for electrical CPR urgent use





CHOICE OF COLOURS

The practical, modest and very attractive design of the Eleganza 1 bed together with a wide range of colour effects creates a pleasant and up-to-date hospital environment.



Safe working load

Castors		Battery	
Dimensions	Braking	Without battery	With battery
125 mm	Individually braked	185 kg	
	Centrally braked	200 kg	250 kg
150 mm	Centrally braked	200 kg	250 kg

Technical parameters

	Castors	
	125 mm	150 mm
External dimension	218 × 99 cm	
Bed extension	15 cm	
Mattress platform	4 section, electric adjustable	
Mattress platform dimensions	200 × 90 × 14 cm	
Electric height adjustment	37 – 73.5 cm	39.5 – 76 cm
Electric backrest adjustment	0 – 70°	
Electric thighrest adjustment	0 – 34°	
TR / ATR	15° – 15°	
Ergoframe®	10 + 6 cm	
Bed weight (as specified)	135 kg	
Standard	IEC: 60601-2-52	
Protection against water and dust	IPx4	
Overload protection	Yes	
Clearance for server	16 cm	18.8 cm



Eleganza 1

Patient

- Ergoframe® against pressure ulcers
- Easy mobilisation
- Available with a number of controls



- Easy clean concept
- Safety features such as the Extender or Protector
- Easy bed manoeuvrability



Staff

Management

- Low TCO – total cost of ownership
- Universal bed for a low price
- Overload protection



 wissner-
bosserhoff

 LINET®

официальный дистрибьютор LINET
в Российской Федерации
ООО "МЕДЭК Инструмент Групп"
119602, Москва, ул. Коштоянца 20/4
тел.: +7(495)755-55-31, 755-55-32
факс +7(495)755-55-30
mail@linet-ru.com, medec@medec.ru
WWW.LINET-RU.COM

Members of LINET Group

LINET spol. s r.o.

ŽelEVčice 5 | 274 01 Slaný | Czech Republic
tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com

